

IMPROVE EFFICIENCIES IN THE OPERATING ROOM AND MINIMIZE TISSUE LOSS


WEISS GLASS CANNULA - LWVI

The tissue has been prepared in an aseptic manner using the modified scuba technique, "F" or "S" stamped (optional), trephinated with a corneal punch, stained with trypan blue and loaded into the glass cannula.

Typical pre-punches are 7.5mm, 7.75mm and 8mm

Approximate Incision Size 2.0-2.2

Glass cannula is packaged ready to attach to a syringe via luer lock connection.

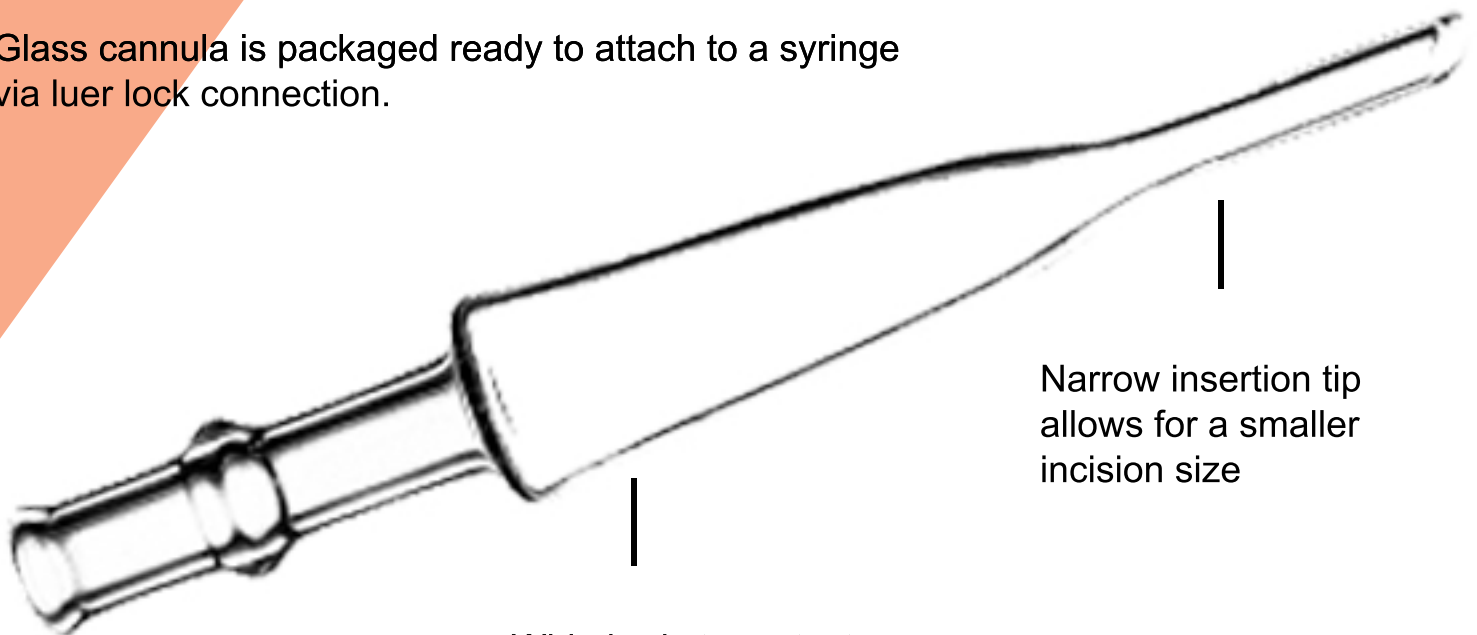

Narrow insertion tip allows for a smaller incision size

Wide body to protect endothelial cells

